



## Installation and Operating Instructions for the iAdapter-Mini™ 2 & 4

Loosen the 8 (eight) screws from the back of the iAdapter™ as shown in the picture below. These screws are captive screws and will not fall out after you loosen them and turn the bottom housing over.

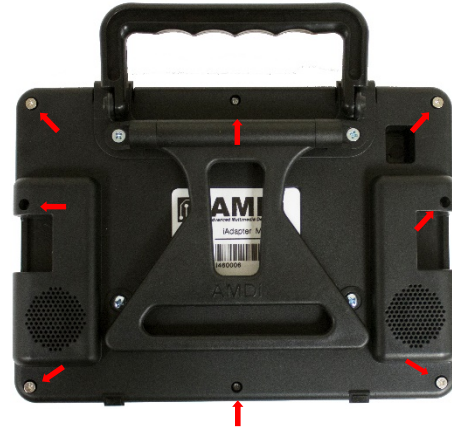

Place the iPad™ onto the 6 (six) rubber supports inside the iAdapter™. The rubber supports are color coded for the different iPad's. RED is for iPad Mini versions 1 through 3. BLACK is for the iPad Mini 4.

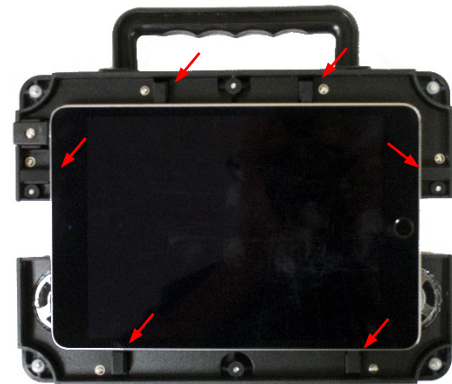

Make sure that the reset button of the iPad® lines up with the reset button of the iAdapter™ as shown in the picture to the right.




Place the top housing of the iAdapter™ over the iPad®. Align the top housing of the iAdapter™ with the bottom housing of the iAdapter™. Carefully hold both the top and the bottom housing of the iAdapter™ together and then turn it over onto a flat surface. Tighten the 8 (eight) screws that will hold the iAdapter™ together.



## Installation and Operating Instructions for the iAdapter-Mini™ 2 & 4

### Charging the battery inside the iAdapter™ and the battery of the iPad™:

The charging jack on the iAdapter™ only charges the battery inside the iAdapter™ for the amplifier. Plug the USB cable supplied with the iAdapter™ into the USB Charger Input jack on the side of the iAdapter™ (see picture at right). Plug the other end of the USB cable into the USB Charger supplied with the iAdapter™ and then plug the charging adapter in a wall socket. The Charger LED on the side of the iAdapter™ will light RED and the LED on the charging adapter will turn GREEN showing power is connected to the charging adapter and the iAdapter. When the Charger LED on the iAdapter changes to GREEN the battery is fully charged. It is recommended that you slide the Power switch to the OFF position while charging the battery of the iAdapter.

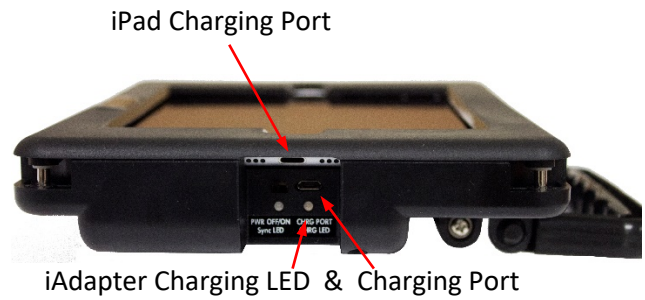

USB Cable

USB Charger